

COLPOSCOPE OP-C2L

THE BEST VISION IN COLPOSCOPY

- LED illumination
- 3 different magnifications
- Light and compact



TECHNOLOGY



producing quality

The OP-C2L colposcope combines the standard qualities of the OPTOMIC colposcopes with the latest LED technology illumination. The LED, integrated in the colposcope head, offers white and uniform lighting, perfect for carrying out the daily work of the specialist.

COLPOSCOPE **OP-C2L**



INCLINED BINOCULAR HEAD




LED ILLUMINATION



OPTIONAL VIDEO SYSTEM

- LED lifespan of 50.000 hours.
- 45° inclined stereoscopic binocular head.
- Optics of great brightness with multilayer treatment.
- Wide field eyepieces, both with dioptic adjustment and protective rubber shells.
- One central handle to operate the head unit easily.
- The arm is an autocompensated parallelogram which can be placed in any position.
- Very soft and smooth operation of all functions.
- High intensity illumination. Built-in green filter in the head.
- 3 steps Galilei magnifications changer and micrometric focusing system.
- Working distance of 250mm with standard objective (f=300mm optional).
- Range of magnifications: 3.5x up to 17.5x, depending on the chosen combination of objective, binocular and eyepieces.
- Different floor and wall mounts available.
- Optional monitor supports and trays for laptop or CCD camera.
- Optional systems for video-colposcopy or digital photography.

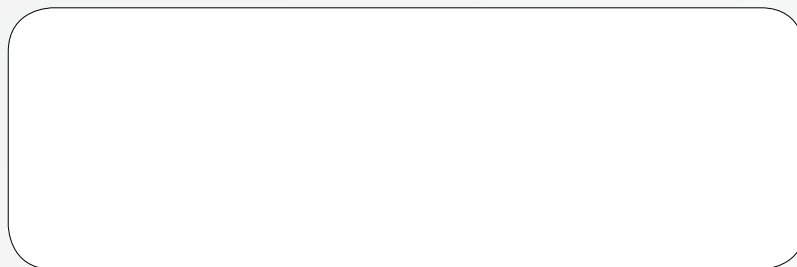

Specifications are subject to change without notice

