





Optimisations at a glance

Overview of all product optimisations made to the medication distribution system

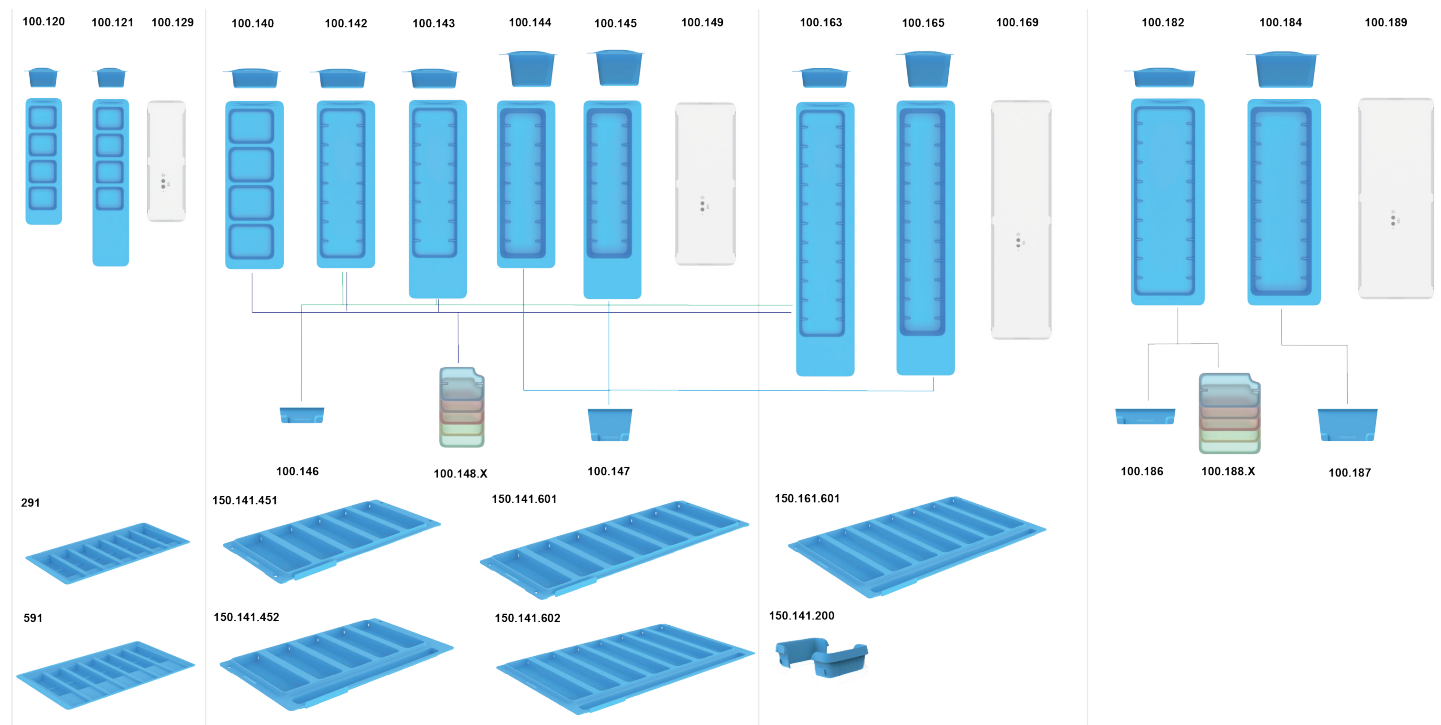

Product range optimisation


The consistent coordination of MediDispensers, dividers and sliding lids has optimised the number of components and the range of variable product dimensions. The diversity of the overall range remains completely intact. In addition, the new medication distribution system enables additional application and combination options.

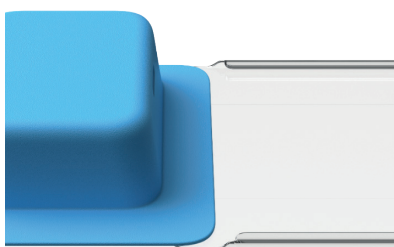

Sliding geometry simplifies the attachment of the sliding lid




The stop bar allows the MediDispenser to be opened in one direction only

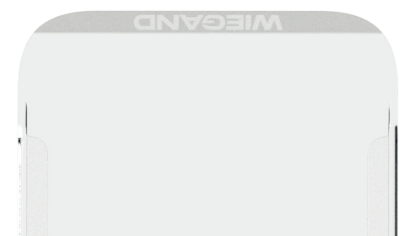


Position line with logo improves visibility of sliding lid

Improved sliding lid

The sliding lid now features a sophisticated sliding profile geometry, which significantly improves the sliding mechanism. The position line with logo, in combination with the dividers, improves the visibility of the sliding lid. The stop bar on the MediDispenser helps to ensure that the medication sequence is followed correctly.

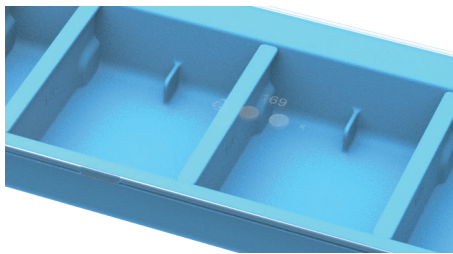

Newly developed plug-in rib for dividers

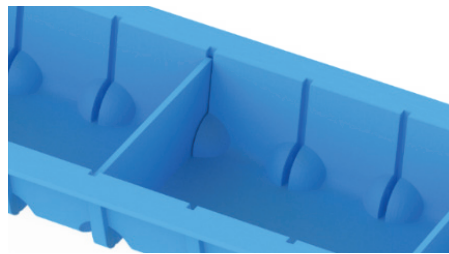


Current plug-in mechanism




Improved grip and a wider top edge for better visibility

Product design

The product design of the MediDispenser has been refined to improve hygiene, and the plug-in ribs for the dividers have been simplified. There are no longer any joints in which dirt can accumulate, and the compartment sizes and spacing of the plug-in ribs have been standardised to 25 mm to improve clarity during the preparation process.

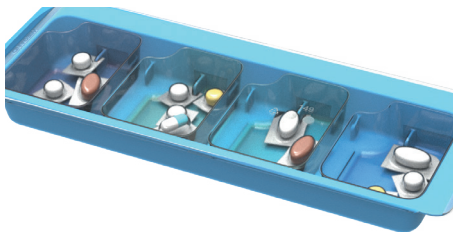

Inlay dishes are available in five colours

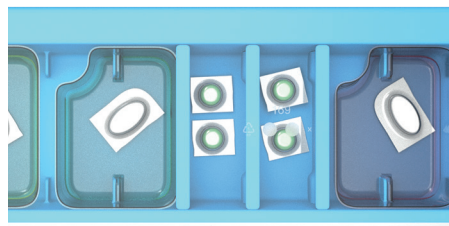

Inlay dishes can be used flexibly




Inlay disher in flat MediDispensers/
MediCups in deep MediDispensers

New combination possibilities

The MediDispensers can be used in a variety of ways and freely combined with distribution aids such as inlay dishes or MediCups. The inlay dishes fit into the flat MediDispensers and can be freely combined with dividers. MediCups can also be used upright in deep MediDispensers.



Fully compatible with WIEGAND FlexLine
trolleys

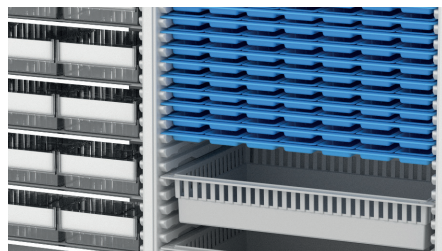

Space-saving storage in the cabinet



Can be used in Wiegand distribution trays

Compatibility with existing Wiegand solutions

The new MediDispensers are compatible with the existing range (some restrictions may apply). Distribution trays can continue to be used without restriction. Most storage and weekly trays remain usable. The new MediDispensers will have slightly more width clearance on old trays. However, the old MediDispensers do not fit on new week and storage trays. Storage in cabinets and transport in WIEGAND FlexLine trolleys with profile walls are still possible without restriction.