




Baby scale "Easy" with fold function

8320

- ▶ Foldable baby scale with integrated carrying handles
- ▶ Lightweight and handy
- ▶ Eight-point measuring system guarantees exact weight evaluation
- ▶ Eight load cells
- ▶ Easy and fast opening: Pull the handles apart and remove the base

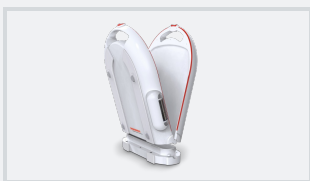

Features

The first baby scale with a handy fold function for parents, midwives and pharmacies can be reduced easily to the half of its length and is impressive because of its low weight. Furthermore the scale is extra handy and easy to transport thanks to the integrated carry handle and a beauty case-shaped design. Even when the child is very fitful the eight-point-weighing guarantees an every time accurate determination of weight which is completed by an easy to read display. The baby scale has an environment-friendly auto-off function for a long battery life.

- ▶ Patented folding mechanism with magnetic lock
- ▶ Low space requirement, easy to transport
- ▶ Easy to read LC display, 15 mm digit height
- ▶ Tare function: Covers/paddings can be tared
- ▶ Hold function for easy reading of the weight value after weighing
- ▶ Switchable to kg/lb units
- ▶ Energy-saving auto-off function
- ▶ Battery operation with two AAA alkaline batteries (included in scope of delivery)

Products

Order number Variation	Approvable	Division g	Maximum load kg
8320.01.001 Baby scale "Easy" with fold function	No	10	15

Accessories

Order number
Variation

Description

5040.01.001



Carry bag black, with flexible shoulder strap and two carrying handles
