


Quality time for better care

"Quality time for better care" is the brand promise of Terumo Medical Care Solutions

HOSPITAL CARE

Surflo™ -W

Winged I.V. Catheter
Product Information

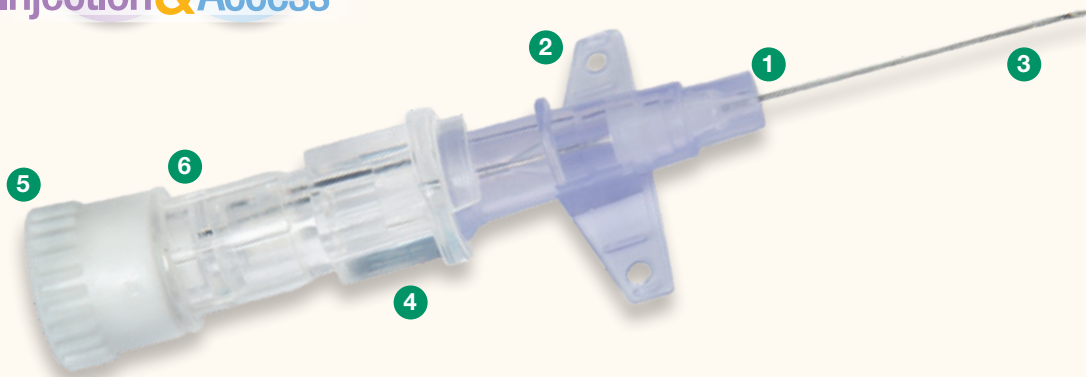

- Ergonomically designed catheter hub – Allows for traditional insertion techniques while providing comfort, ease of stabilization and reduced movement in the vessel
- The flashback chamber visually confirms that the catheter is correctly inserted in the vein
- The sharp needle tip minimizes trauma to the vessel during use
- Full range from 14G to 26G

Description	Raw Materials
1 Catheter	FEP + Ba SO ₄
2 Catheter hub Caulking pin	Polypropylene + master batch Polyoxymethylene
3 Cannula (internal needle)	Stainless steel
4 Needle hub	Styrene butadiene copolymer
5 Luer-lock cap	Polypropylene
6 Vent fitting filter Vent fitting body	Polyethylene Styrene butadiene copolymer
Protector	Polyethylene
Unit pack	Adhesive coated lid paper. Thermoformed blister film: laminated polyethylene/ polypropylene complex

Does not contain components made of PVC or natural rubber latex.

Designed and
manufactured by skilled
technicians in Europe

Your catheters of choice for a wide range of I.V. applications

Product Range

Product Code	Colour	External Catheter Ø		Internal Catheter Ø	Catheter Length	Flow Rate ml/min.
SR+DM2619WX	Violet	26G	0.64 mm	0.45 mm	19 mm	19
SR+DM2419WX	Yellow	24G	0.74 mm	0.55 mm	19 mm	29
SR+DM2225WX	Blue	22G	0.90 mm	0.65 mm	25 mm	42
SR+DM2032WX	Pink	20G	1.00 mm	0.75 mm	32 mm	59
SR+DM1832WX	Green	18G	1.30 mm	0.95 mm	32 mm	103
SR+DM1845WX	Green	18G	1.30 mm	0.95 mm	45 mm	96
SR+DM1745WX	White	17G	1.50 mm	1.15 mm	45 mm	155
SR+DM1645WX	Grey	16G	1.75 mm	1.35 mm	45 mm	225
SR+DM1445WX	Orange	14G	2.00 mm	1.55 mm	45 mm	290



Packaging

Single Unit Blister Pack

Unit Box 50 pieces

Shipping Carton 500 pieces, Gross Weight ± 3.5 kg
Dimensions 390 x 730 x 140 mm

Regulatory Status

CE Certification – Class IIa medical device

Notified body
TÜV Rheinland (0197)

Manufacturer
TERUMO EUROPE N.V., Interleuvenlaan 40, 3001 Leuven, Belgium

Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10⁻⁶

Storage

- Shelf life of 5 years
- The batch number and expiry date are printed on the paper of the unit pack
- Do not store at extreme temperature and humidity

Precautions

- Do not use if unit package is damaged
- Do not re-use / do not re-sterilize
- Do not attempt to re-insert a partially or completely withdrawn needle
- Dispose of safely after single use to avoid risk of infection
- Application time should be consistent with accepted standards of practice

Always consult the product label for a complete overview of warnings, cautions and/or precautions prior to actual use.

This range complements the wider Terumo product offering including a comprehensive portfolio of both safety and conventional I.V. products which can be found at our web page on www.terumo-europe.com

Your Terumo partner



TERUMO EUROPE N.V.

Terumo Corporation Terumo Europe N.V.
+81 3 6742 8500 +32 16 38 12 11

www.terumo-europe.com

All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners. Published by Terumo Europe N.V.