


**HOSPITAL CARE**

# Surflo™ S-PLUS

Straight I.V. Catheter  
Product Information

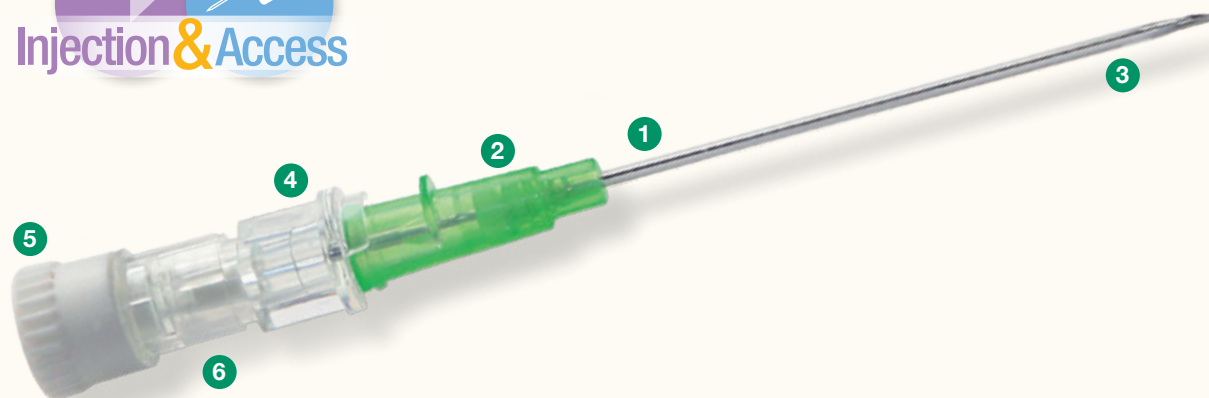

- Ergonomically designed catheter hub – Allows for traditional insertion techniques while providing comfort, ease of stabilization and reduced movement in the vessel
- The flashback chamber visually confirms that the catheter is correctly inserted in the vein
- The sharp needle tip minimizes trauma to the vessel during use
- Full range from 14G to 26G

Description	Raw Materials
<b>1</b> Catheter	FEP + Ba SO <sub>4</sub>
<b>2</b> Catheter hub Caulking pin	Polypropylene + master batch Polyoxymethylene
<b>3</b> Cannula (internal needle)	Stainless steel
<b>4</b> Needle hub	Styrene butadiene copolymer
<b>5</b> Luer-lock cap	Polypropylene
<b>6</b> Vent fitting filter Vent fitting body	Polyethylene Styrene butadiene copolymer
Protector	Polyethylene
Unit pack	Adhesive coated lid paper. Thermoformed blister film: laminated polyethylene/ polypropylene complex

Does not contain components made of PVC or natural rubber latex.

Designed and  
manufactured by skilled  
technicians in Europe

## Your catheters of choice for a wide range of I.V. applications

### Product Range

Product Code	Colour	External Catheter Ø		Internal Catheter Ø	Catheter Length	Flow Rate ml/min.
SR+DM2619SX	Violet	26G	0.64 mm	0.45 mm	19 mm	19
SR+DM2419SX	Yellow	24G	0.74 mm	0.55 mm	19 mm	29
SR+DM2225SX	Blue	22G	0.90 mm	0.65 mm	25 mm	42
SR+DM2032SX	Pink	20G	1.00 mm	0.75 mm	32 mm	59
SR+DM1832SX	Green	18G	1.30 mm	0.95 mm	32 mm	103
SR+DM1845SX	Green	18G	1.30 mm	0.95 mm	45 mm	96
SR+DM1745SX	White	17G	1.50 mm	1.15 mm	45 mm	155
SR+DM1645SX	Grey	16G	1.75 mm	1.35 mm	45 mm	225
SR+DM1445SX	Orange	14G	2.00 mm	1.55 mm	45 mm	290



### Packaging

Single Unit Blister Pack

Unit Box 50 pieces

Shipping Carton 500 pieces, Gross Weight ± 3.0 kg  
Dimensions 390 x 730 x 140 mm

### Regulatory Status

CE Certification – Class IIa medical device

Notified body  
TÜV Rheinland (0197)

Manufacturer  
TERUMO EUROPE N.V., Interleuvenlaan 40, 3001 Leuven, Belgium

### Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10<sup>-6</sup>

### Storage

- Shelf life of 5 years
- The batch number and expiry date are printed on the paper of the unit pack
- Do not store at extreme temperature and humidity

### Precautions

- Do not use if unit package is damaged
- Do not re-use / do not re-sterilize
- Do not attempt to re-insert a partially or completely withdrawn needle
- Dispose of safely after single use to avoid risk of infection
- Application time should be consistent with accepted standards of practice

Always consult the product label for a complete overview of warnings, cautions and/or precautions prior to actual use.

This range complements the wider Terumo product offering including a comprehensive portfolio of both safety and conventional I.V. products which can be found at our web page on [www.terumo-europe.com](http://www.terumo-europe.com)

Your Terumo partner



#### TERUMO EUROPE N.V.

Terumo Corporation Terumo Europe N.V.  
+81 3 6742 8500 +32 16 38 12 11

[www.terumo-europe.com](http://www.terumo-europe.com)

All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners. Published by Terumo Europe N.V.