


HOSPITAL CARE

Surshield™ Versatus -WP

Winged and Ported I.V. Catheter with Protector
Product Information

The Surshield protector allows to completely close the needle, in order to prevent accidental needle stick injuries when disposing of the intravenous cannula. Using a safety catheter does not require any change in the insertion technique.

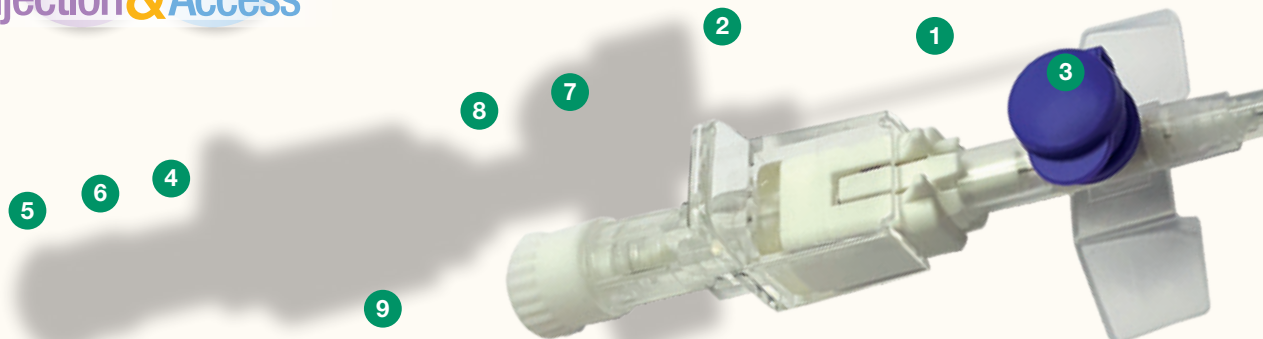

| Description | | Raw Materials |
|----------------------------|------------------------------------------|-------------------------------------------------|
| 1 | Catheter | Polyurethane + BaSO ₄ |
| 2 | Winged Catheter hub Caulking pin | Polypropylene Polyoxymethylene |
| 3 | Cannula (internal needle) | Stainless steel |
| 4 | Needle hub | Acrylonitrile butadiene styrene |
| 5 | Luer-lock cap | Polypropylene |
| 6 | Vent fitting filter Vent fitting body | Polyethylene Acrylonitrile butadiene styrene |
| 7 | Port Closure | Polypropylene + master batch |
| 8 | Injection port valve | Silicone tube |
| | Protector | Polyethylene |
| Surshield Protector | | |
| 9 | Protection housing | Polyamide |
| | Safety fingers | Polysulfone |
| | Sliding bush + Needle blocker | Polycarbonate + stainless steel cup |

Does not contain components made of PVC or natural rubber latex.

Designed and
 manufactured by skilled
 technicians in Europe

Your catheters of choice, part of a wider I.V. portfolio

Product Range

| Product Code | Colour | External Catheter Ø | | Internal Catheter Ø | Catheter Length | Flow Rate ml/min. |
|--------------|--------|---------------------|---------|---------------------|-----------------|-------------------|
| SR+DS2619PD | Violet | 26G | 0.64 mm | 0.45 mm | 19 mm | 19 |
| SR+DS2419PD | Yellow | 24G | 0.74 mm | 0.55 mm | 19 mm | 29 |
| SR+DS2225PD | Blue | 22G | 0.90 mm | 0.65 mm | 25 mm | 42 |
| SR+DS2032PD | Pink | 20G | 1.00 mm | 0.75 mm | 32 mm | 59 |
| SR+DS1832PD | Green | 18G | 1.30 mm | 0.95 mm | 32 mm | 103 |
| SR+DS1845PD | Green | 18G | 1.30 mm | 0.95 mm | 45 mm | 96 |
| SR+DS1745PD | White | 17G | 1.50 mm | 1.15 mm | 45 mm | 155 |
| SR+DS1645PD | Grey | 16G | 1.75 mm | 1.35 mm | 45 mm | 225 |
| SR+DS1445PD | Orange | 14G | 2.00 mm | 1.55 mm | 45 mm | 290 |



Packaging

Single Unit Blister pack

Unit Box 50 pieces

Shipping Carton 200 pieces, Dimensions 390 x 384 x 180 mm

Regulatory Status

CE Certification – Class IIa medical device

Notified body

Istituto Superiore di Sanità 0373

Manufacturer

Delta Med S.p.A. (Italy)

Distributor

TERUMO EUROPE N.V., Interleuvenlaan 40, 3001 Leuven, Belgium

Sterilization

With ethylene oxide to a Sterility Assurance Level of at least 10⁻⁶

Storage

- Shelf life of 5 years
- The batch number and expiry date are printed on the paper of the unit pack
- Do not store at extreme temperature and humidity

Precautions

- Do not use if unit package is damaged
- Do not re-use / do not re-sterilize
- Do not attempt to re-insert a partially or completely withdrawn needle
- Dispose of safely after single use to avoid risk of infection
- Application time should be consistent with accepted standards of practice

Always consult the product label for a complete overview of warnings, cautions and/or precautions prior to actual use.

This range complements the wider Terumo product offering including a comprehensive portfolio of both safety and conventional I.V. products which can be found at our web page on www.terumo-europe.com

Your Terumo partner



TERUMO EUROPE N.V.

Terumo Corporation Terumo Europe N.V.
+81 3 6742 8500 +32 16 38 12 11

www.terumo-europe.com

All brand names are trademarks or registered trademarks of TERUMO CORPORATION and their respective owners. Published by Terumo Europe N.V.