



producing quality

EXAMINATION CHAIR

OP-MH SERIES

**OPTIMAL DESIGN FOR
GYNECOLOGY AND OBSTETRICS**

- Exceptional minimum height:
down to 48cm
- High load capacity: 310 kg
- Motorized leg rest



OP-MH A CHAIR DESIGNED FOR EFFICIENCY

The OP-MH and OP-MH+ chairs with extendable leg rest are designed to meet all gynecology and obstetrics requirements.

Their versatility comes from a wide range of movements powered by four independent motors, along with a modular accessory system that allows professionals tailor each chair to the specific needs of their practice.

Manufactured with the highest quality materials, paying close attention to ergonomic needs, and using a comfortable and hygienic padding system, OPTOMIC ensures the best experience for both doctor and patient.




BARRIER-FREE ACCESS FOR PATIENTS AND PROFESSIONALS

OPTOMIC's expertise is reflected in the OP-MH chairs, designed for optimal specialist efficiency and patient comfort.

■ MINIMUM HEIGHT OF UP TO 48cm*

Engineered for effortless patient access, with special attention to elderly patients, patients with reduced mobility, and those with disabilities.

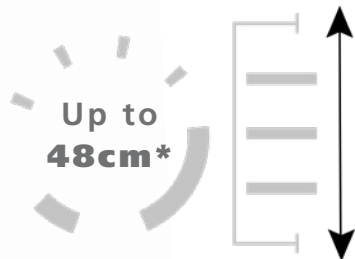

* Only on OP-MH+ model



■ LEG REST MOVEMENT FOR PATIENT CARE

The motorized leg rest enables the patient to reach the correct examination position, without assistance from the healthcare professional.



■ LIFTING CAPACITY UP TO 310 kg

Professionals can confidently use the chair with heavier patients, knowing it meets the most demanding requirements.

■ EASY-ACCESS FOOT SUPPORTS

The low minimum height facilitates easy patient entry and exit.



INDEPENDENT AND MOTORIZED MOVEMENTS

Easily set the desired position and save it in one of the four internal memories for the next examination. The independent movement of the OP-MH elements allows it to adapt to the specific needs of each specialty.

TOTAL CONTROL THROUGH HAND CONTROL OR FOOT PEDAL

All movements can be controlled using these two devices.

- Hand control: All movements, direct access to Trendelenburg, lock button, 3 programmable memories, and 4 preprogrammed memories.
- Foot pedal: All movements and 3 programmable memories.



MOTORIZED RETRACTABLE LEG REST

The leg rest folds under the seat, providing the professional with ample space to move and work comfortably.

TRENDELENBURG POSITION



A shortcut button on the hand control quickly positions the chair in the Trendelenburg position when required.



"GYNECOLOGY" POSITION

GYN

This pre-programmed memory automatically positions the chair for optimal gynecological examinations.



"OBSTETRICS" POSITION

OBS

This pre-programmed memory sets the chair to the stretcher position, ensuring a comfortable obstetric evaluation.



EXAMINATION CHAIR **OP-MH** series

INDEPENDENT BACK AND LEG REST MOVEMENT

Complete control of all movements through the hand control or foot pedal.

The combined adjustability of backrest, seat, and leg rest allows the OP-MH chair to meet all postural requirements in gynecology and obstetrics.

PELVIC MOVEMENT

Facilitates examinations by allowing the chair to be placed in a comfortable and natural position both for the professional and the patient.




+22°

-6°



UP/DOWN MOVEMENT

Exceptionally low minimum height to facilitate patient access, especially for those with reduced mobility.

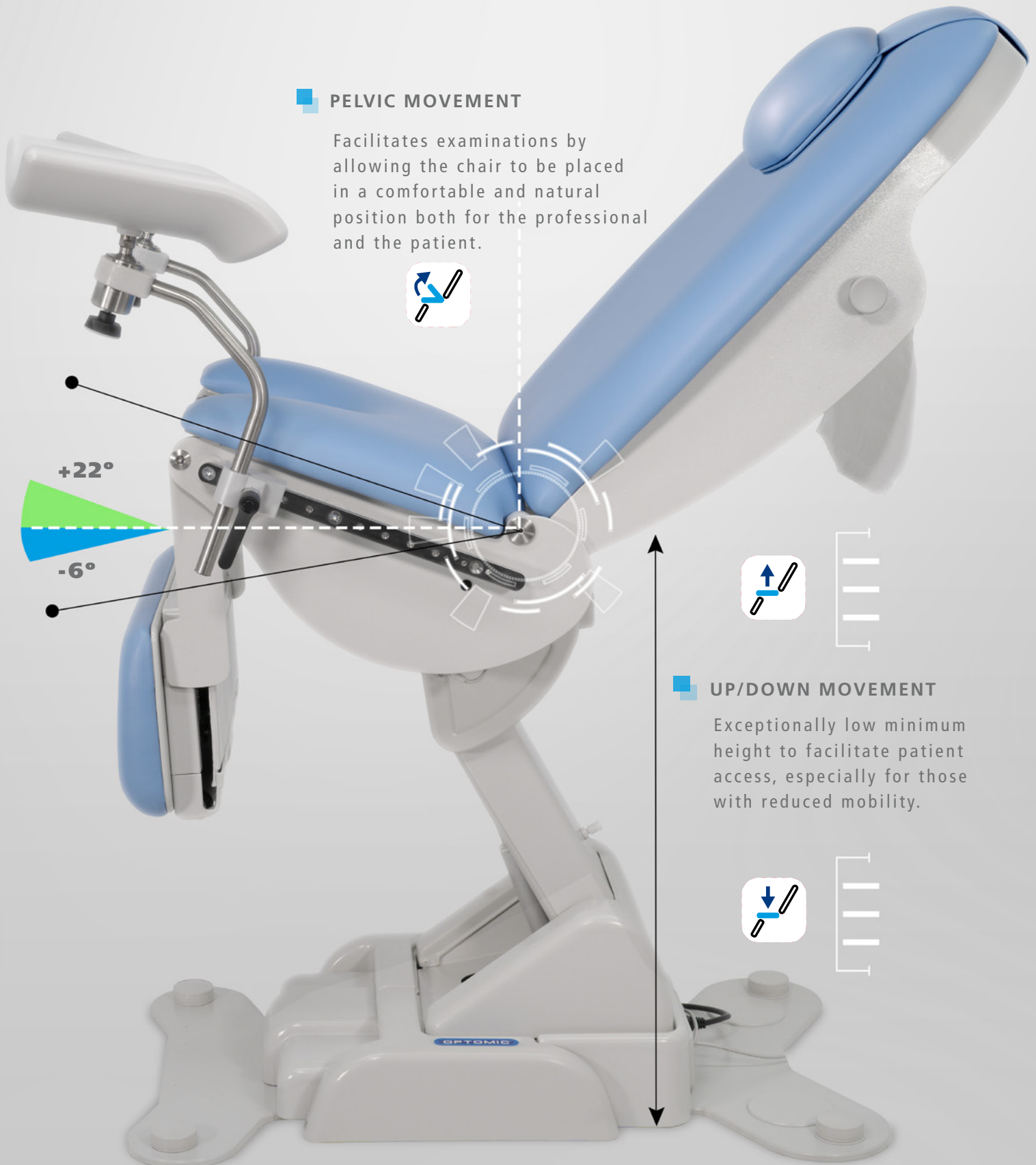

EXAMINATION CHAIR **OP-MH** series

ACCESSORIES FOR EACH AND EVERY NEED

Thanks to its side rail system, the OP-MH can be easily configured with any of our accessories to meet the specialist's requirements. Footrests, leg rests, and armrests can be adjusted to ensure maximum comfort for both the professional and the patient.

LEG SUPPORTS WITH FOOT SUPPORTS

A 3-position system that adjusts the foot support angle for optimal ergonomics. It also allows lateral opening with a gradual brake system.

Magnetically fixed cushioned support for barefoot patients. The support is easily removed for patients wearing shoes.

GOEPEL LEG SUPPORTS

Made of stainless steel and padded with injected flexible polyurethane, very durable and easily washable.

The OP-MH leg supports can be adjusted to any position in order to reach the best examination position.



Foot support option:
Adjustable




Foot support option:
Fixed



Goepel option:
360° rotatable with
complete locking



Goepel option: With
adjustable pan-tilt
mechanism

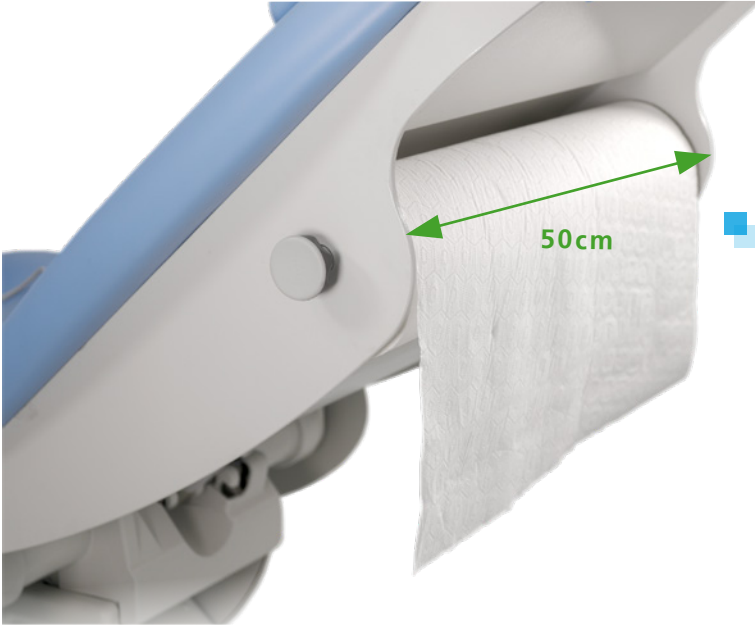


PAPER ROLL HOLDER

The chair includes a roll holder to place the paper on seat and backrest.

The integrated design saves space and simplifies access and paper roll replacement. It accommodates standard 50 cm rolls, and, with the optional extension, it can hold 60 cm rolls

AUXILIARY CUSHION (Optional)

Padded cushion with adjustable height and magnetic placement

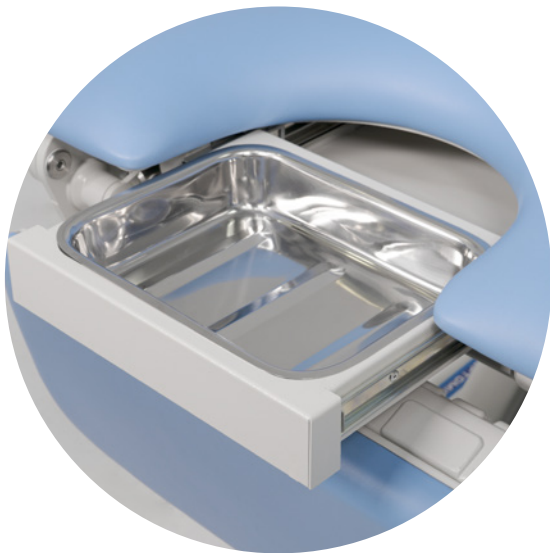

SECRETION TRAY

Made of stainless steel and easy to remove and clean, 2.5 L capacity



COLPOSCOPE SUPPORT (Optional)

OP-MH chairs can be fitted with a wide range of colposcopes using a sturdy support specifically designed for this purpose.



ARMRESTS
(Optional)



PAPER ROLL HOLDER
EXTENSION UP TO 60cm
(Optional)



HIDDEN WHEEL SYSTEM
(Optional)

FEATURES

MOVEMENT:

LOWERING/LIFTING

Motorized column: OP-MH+: 48 cm min.- 102 cm max.

OP-MH: 55 cm min.- 102 cm max.

BACKREST

Through electric motor. -22° min. / 70° max.

LEG REST

Through electric motor. 0° min. / 110° max.

PELVIC

Through electric motor. 20° max. above/-6° min. below horizontal

TRENDELENBURG

Tilting -20°

ANTI-TRENDELENBURG

Tilting +6°

HAND CONTROL AND FOOT PEDAL

Control of all movements.

MEMORIES

Seven memory settings: 3 programmable and 4 pre-programmed.

LEG SUPPORTS

Two types of Goepel (Standard and Ergo)

Two types of Foot Support (Standard and Ergo)

SECRETIONS TRAY

Made of stainless steel and removable. 2,5L capacity

ROLL HOLDER

Paper reel 50cm max.

UPHOLSTERY

Washable and antibacterial PUR-PVC Polyester. Wide range of colors.

FIREPROOF STANDARD

Flame retardant class M2 - BS5852 1-2

MAINS CONNECTION

110 /240V / 50-60Hz

MAXIMUM LOAD

310 kg

OPTIONAL ACCESSORIES

COLPOSCOPE SUPPORT

Compatible with the main brands on the market.

Colposcope can be placed at the right or left side of the chair.

ARMREST

Removable. Adjustable in height and longitudinal position.

ROLL HOLDER EXTENSION

Accommodates paper rolls up to 60cm.

AUXILIARY CUSHION

Padded cushion with adjustable height and magnetic placement

EXAMINATION LAMP

High power LED with dimmer switch and flexible arm.

WHEEL SYSTEM

Hidden wheel system that can be unfolded by means of a lever.

OPTOMIC ESPAÑA S.A.

Madroño 4. E-28770, Colmenar Viejo. Madrid - Spain

Tel: +34 91 804 44 44

export@optomic.com

www.optomic.com



producing quality



Specifications are subject to change without prior notice

